




Stille osteotomes, gouges and mallets

Stille osteotomes, gouges and mallets

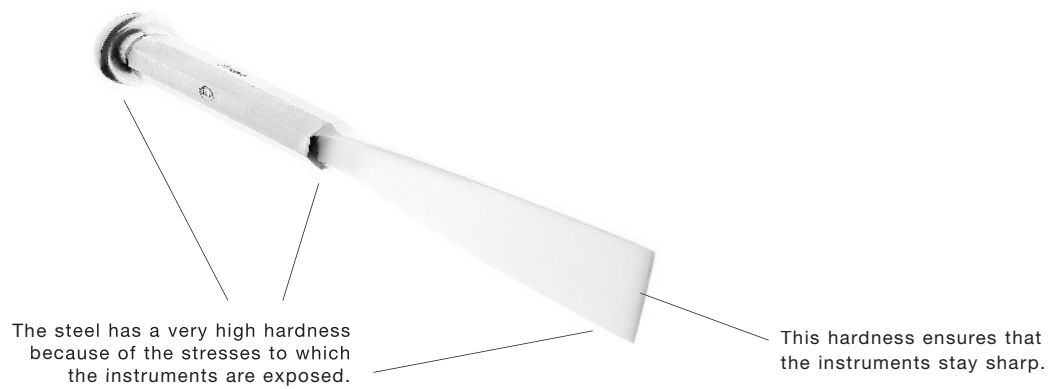

Stille osteotomes, gouges and mallets are made of choice quality steel with high HRC hardness. Stille was the first company in the world to manufacture osteotomes and we are today world famous for the superior cutting performance. The special grinding of Stille osteotomes results in superior quality instruments with very strong tips that are prevented from getting stuck in the bone. The phenolic and stainless steel handles are designed to prevent from slipping.


Osteotomes and gouges must be very hard to withstand the stress they are exposed to. At the same time it is important that the material gives way when touching implants so that these will not be damaged. The hardness also ensures that the instruments stay sharp for a long period of time.



STILLE-KAYE

Delicate osteotome with cross-serrated head and comfortable handle. Available in 7 different sizes.


* 111-21306-2	2 mm	16.5 cm 6.5"
111-21306-4	4 mm	16.5 cm 6.5"
111-21306-6	6 mm	16.5 cm 6.5"
111-21306-8	8 mm	16.5 cm 6.5"
111-21306-10	10 mm	16.5 cm 6.5"
111-21306-12	12 mm	16.5 cm 6.5"
111-21306-15	15 mm	16.5 cm 6.5"

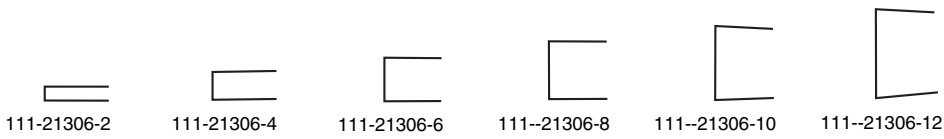

MINI LAMBOTTE

Delicate mini osteotome with rectangular handle. 8 different sizes.

111-21450-2	Straight	2 mm	12.5 cm 5"
111-21450-3	Straight	3 mm	12.5 cm 5"
111-21450-4	Straight	4 mm	12.5 cm 5"
111-21450-5	Straight	5 mm	12.5 cm 5"
* 111-21450-6	Straight	6 mm	12.5 cm 5"
111-21450-8	Straight	8 mm	12.5 cm 5"
111-21450-10	Straight	10 mm	12.5 cm 5"
111-21450-12	Straight	12 mm	12.5 cm 5"

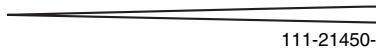

STILLE OSTEOTOME

Delicate osteotome with a comfortable phenolic handle.

* 111-21456-2	Straight	2 mm	16 cm 6.25"
111-21456-4	Straight	4 mm	16 cm 6.25"
111-21456-6	Straight	6 mm	16 cm 6.25"
111-21456-8	Straight	8 mm	16 cm 6.25"
111-21456-10	Straight	10 mm	16 cm 6.25"
111-21456-12	Straight	12 mm	16 cm 6.25"



OSTEOTOMES


* Photo shows this instrument. Not shown at actual size



CINELLI

Double-guarded osteotome for hump removal.




- * 111-21516-14 14 mm 16 cm 6.25"
- 111-21516-16 16 mm 16 cm 6.25"

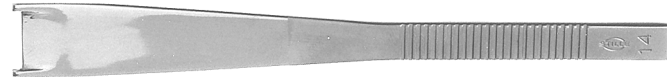

SILVER OSTEOTOME

Single-guarded osteotome for lateral osteotomy.

- 111-21500-1 Straight 18 cm 7"
- * 111-21500-2 Curved right 18 cm 7"
- 111-21500-3 Curved left 18 cm 7"

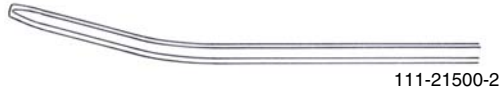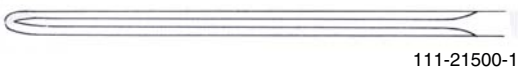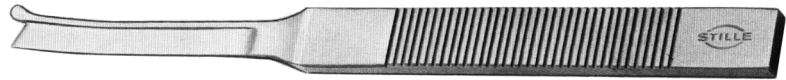

NEIVERT OSTEOTOME

Delicate single-guarded osteotome for lateral osteotomy.

- 111-21501-1 Curved left 19.5 cm 7.75"
- * 111-21501-2 Curved right 19.5 cm 7.75"
- 111-21501-3 Straight 19.5 cm 7.75"



OSTEOTOMES

* Photo shows this instrument. Not shown at actual size.




STILLE ORIGINAL GOUGE

Original Stille gouge. Available in five different sizes from 6 mm to 15 mm.



* 111-21315-6	6 mm	16.5 cm	6.5"
111-21315-8	8 mm	16.5 cm	6.5"
111-21315-10	10 mm	16.5 cm	6.5"
111-21315-12	12 mm	16.5 cm	6.5"
111-21315-15	15 mm	16.5 cm	6.5"

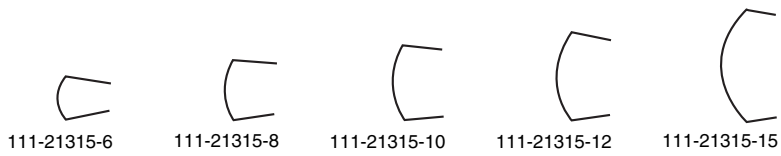

STILLE MALLET

Stainless steel mallet.

* 111-21820-	200 g	7 oz	21.5 cm	8.5"
--------------	-------	------	---------	------



STILLE MALLET WITH PHENOLIC HANDLE

Mallet with stainless steel head, phenolic handle and interchangeable teflon insert.

* 111-21825-	145 g	5 oz	17 cm	6.75"
111-21825-5	Teflon spare insert			



G
O
U
G
E
S
&
M
A
L
L
E
T
S


* Photo shows this instrument. Not shown at actual size



STILLE MALLET

Mallet with plastic head and stainless steel handle.

- * 111-21920- 245 g 8.5 oz diam. 40 mm 23.5 cm 9.25"
- 111-21920-5 Spare plastic head

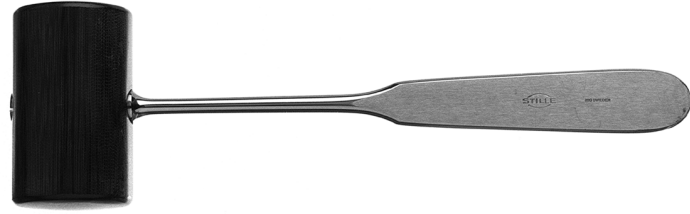

M
A
L
L
E
T
S

* Photo shows this instrument. Not shown at actual size.